

英文成册, 材质:封面光面128g铜版纸, 内页普通80g书写 纸, 尺寸:100×140mm

# Blood Pressure Monitor

Model Number: F1701T

## USER'S MANUAL



Shenzhen Jamr Technology Co., Ltd.  
A101-301, D101-201, Jamr Science & Technology Park, No. 2  
Guiyuan Road, Guixiang Community, Guanlan Street, Longhua  
District, 518100 Shenzhen, PEOPLE'S REPUBLIC OF CHINA



Shanghai International Holding Corp. GmbH (Europe)  
Eiffestrasse 80, 20537 Hamburg, Germany

Published Date:2022/5/27

Version Number:1.0



**CONTENTS**

1 Introduction and Intended Use.....2

2 Important Information on Blood Pressure and its Measurement.....5

3 Components of Your Blood Pressure Monitor.....6

4 Using Your Monitor for the First Time.....8

5 Measurement Procedure.....10

6 Care and Maintenance.....17

7 Warranty/Service.....18

8 Certifications.....18

9 Technical Specifications.....19

10 EMC Declaration.....20

## 1 Introduction and Intended Use

The device is a fully automatic digital blood pressure measuring device using oscillometric technique to measure systolic and diastolic blood pressure as well as the pulse for adults that ages are more than 12 years old by wrapping around the upper arm with cuff circumference ranging from 22 cm to 42 cm. The device can be used in medical facilities or at home, and only for indoor use.

**Contraindication:** The device is not used for patients under dialysis therapy or on anticoagulant, antiplatelets, or steroids.

The device is provide accurate blood pressure measurement values that are effective and suitable for clinical and home use.

Before using, please read this instruction manual carefully and then keep it in a safe place.

### 1.1 Remember...

- Only a health-care professional is qualified to interpret blood pressure measurements.
- This device is NOT intended to replace regular medical checkups.
- The device is intended for use by adults only and not intended for use on children and pregnant patient. The effectiveness has not been established in pregnant (including pre-eclamptic) patients.
- In cases of irregular heartbeat, measurements made with this instrument should only be evaluated after consultation with a physician.
- The products, including accessories, shall be processed in accordance with local regulations after reaching the life cycle.
- Any serious incident that has occurred in relation to the device should be reported to the manufacturer and the competent authority of the Member State in which the user and/or patient is established.

## 1.2 Warnings and Precautions

**Warning:** Do not use the AC adapter if the unit or the power cord is damaged. Turn off the power and unplug the power cord immediately.

**Warning:** The device is not suitable for use in the presence of flammable anesthetic mixtures with air or with oxygen or nitrous oxide.

**Warning:** If the patient is an intended operator, the functions of monitoring blood pressure and pulse rate can be safely used by patient. The routine clean and changing batteries can be performed by the patient.

**Warning:** The device provides a DC input port connected to external ac adapter. It is recommended that use the adapter specified by the manufacturer. The adapter should meet the following conditions: class II equipment, output voltage: DC 5V, current:  $\geq 1A$ , and comply with IEC 60950, IEC 60601-1 or IEC 62368-1, provide at least two MOOP insulation between ac input and dc output. External adapter connected to medical electrical equipment through the DC input port must comply with the respective IEC or ISO standards (e.g. IEC 60950 or IEC 62368-1 for data processing equipment). Further more all configurations shall comply with the requirements for medical electrical systems (see IEC 60601-1-1 or clause 16 of the 3Ed. of IEC 60601-1, 60601-1-2, respectively). Anybody connecting external adapter to medical electrical equipment configurations a medical system and is therefore responsible that the system complies with the requirements for medical electrical systems. Attention is drawn to the fact that local laws take priority over the above mentioned requirements. If in doubt, consult your local representative or the technical service department.

**Warning:** Too frequent measurements can cause injury to the PATIENT due to blood flow interference.

**Warning:** Don't place the cuff over wound part.

**Warning:** Pressurization of the CUFF can temporarily cause loss of function of simultaneously used monitoring ME EQUIPMENT on the same limb.

**Warning:** Regularly checking the operation of the blood pressure monitor to ensure that it does not cause long-term damage to the patient's blood circulation.

Warning: Apply CUFF and its pressurization on the side of the patient's mastectomy or lymph node removal can cause injury.

Warning: To avoid any possibility of accidental strangulation, keep this device away from children and do not drape tubing around your neck.

Caution: The user must check that the equipment functions safely and see that it is in proper working condition before being used.

Caution: To avoid damaging the device, keep this unit away from children and pets.

Caution: The standard material used for the bladder and tubing is latex-free.

Caution: The device is intended to monitor, not to diagnose. Unusual values have to be always discussed with a physician. Under any circumstance, you should not alter the dosages of any drugs prescribed by a physician.

Caution: The device cannot be used to substitute the professional ECG monitor device for monitoring the frequency of heart beat.

Caution: This device can not be used together with HF surgical equipment.

Note: To obtain the greatest accuracy from your blood pressure monitor, it is recommended that the instrument be used within the specified temperature and the relative humidity, please see the Technical Specifications.

Note: The device can not be used in MRI environment.

Note: The cuff is defined as the applied part. The user should contact the manufacturer for assistance, if needed, replace, or maintaining the device.

Note: This device contains sensitive electronic components. Avoid strong electrical or electromagnetic fields in the direct vicinity of the device (e.g. mobile telephones, microwave ovens) during use. These can lead to erratic results.

Note: Do not attempt to service or repair this device yourself. Should a malfunction occur, refer to local distributor or the manufacturer.

## 2 Important Information on Blood Pressure and its Measurement


### 2.1 How does high or low blood pressure arise?

Your level of blood pressure is determined in the circulatory center of the brain and adjusts to a variety of situations through feedback from the nervous system. To adjust blood pressure, the strength and speed of the heart (Pulse), as well as the width of circulatory blood vessels is altered. Blood vessel width is controlled by fine muscles in the blood vessel walls.







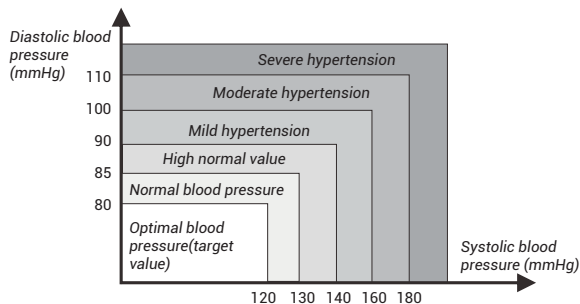































Your level of arterial blood pressure changes periodically during heart activity: During the "blood ejection" (Systole) the value is highest (systolic blood pressure value). At the end of the heart's "rest period" (Diastole) pressure is lowest (diastolic blood pressure value).

### 2.2 Which values are normal?

Please refer to the diagram below (Picture-01)



Picture-01


	WHO grids in device	WHO Classification
                                     		

### 3 Components of Your Blood Pressure Monitor

Diagram illustrating the front panel components of the device:

- Cuff Connector Port**: Located on the left side of the device.
- AC Adapter Port**: Located on the right side of the device.
- LCD Display**: The central screen area.
- SET button**: A small button located below the LCD display.
- Users Switch button**: A button located below the LCD display.
- Memory button**: A button located below the LCD display.
- START/STOP button**: A circular button located below the LCD display.

6

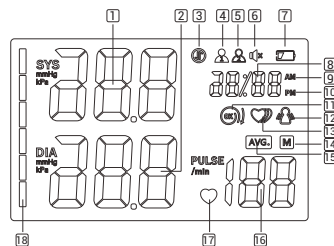

### 3.2 The symbols on the LCD display

- |   |                                  |
|---|----------------------------------|
| 1.Systolic blood pressure                           | 2.Diastolic blood pressure       |
| 3.Bluetooth symbol                                  | 4.USER 1                         |
| 5.USER 2  | 6.Mute symbol                    |
| 7.Low battery symbol                                | 8.Date/Time display              |
| 9.AM symbol   | 10.PM symbol                     |
| 11.Cuff self-checking function                      | 12.Movement error symbol         |
| 13.Irregular heartbeat symbol                       | 14.Memory symbol                 |
| 15.Average value symbol                             | 16.Pulse display / Memory number |
| 17.Heartbeat symbol<br>(Flashes during measurement) | 18.WHO symbol                    |

### 3.3 Features of Model F1701T

1. Voice function
2. Double users: 2 x 120 sets memory
3. Cuff self-checking function
4. Irregular heartbeat checking
5. Average value function
6. Low battery display
7. WHO function
8. Auto power-off
9. External power adapter support
10. Volume adjustment
11. Date/time display
12. Bluetooth function

## 4 Using your Monitor for the First Time


### 4.1 Activating the pre-installed batteries

#### Battery Installation

Use only 1.5V "AA" alkaline batteries with this device.

1. Press the hook on the bottom of the battery cover and lift the cover off in the direction of the arrow (Picture-04).

2. Install 3 "AA" size batteries and make sure the + (positive) and - (negative) polarities match the polarities of the battery compartment, then close the battery cover. Make sure that the battery cover is securely in position.

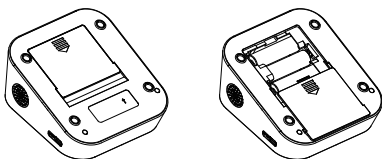

Picture-04

#### Battery replacement

##### Low Battery Indicator

1. When the Low Battery Indicator appears on the display, turn the monitor off and remove all the batteries. Replace with 3 new batteries at the same time. Long-life alkaline batteries are recommended.

2. To prevent the damage of monitor from leaked battery fluid, please take out of battery if the monitor unused in a long time (generally more than 3 months). If battery fluid should get in your eyes, immediately rinse with plenty of clean water. Contact a physician immediately.

3. Dispose of the device, components and optional accessories according to applicable local regulations. Unlawful disposal may cause environmental pollution.

### 4.2. System Settings

Before setting, ensure that all batteries are loaded correctly.

#### Setting the User ID(1 or 2):

With the unit off, Press the (1↓) button and then you can set the User ID user by pressing the (8) button.

#### Setting the Year/ Month/ Date/ Hour system/ Time/ Volume

With the unit off, Long press the (1↓) button for more than 3s, and then you can start to set.

##### Setting the Year

When the year display is flashing, press the (8) button continuously and it will increase continuously 1 by 1 until 2049, and then return the original year, press the (8) button continuously and it will reduce continuously 1 by 1, once the year set is OK, press the (1↓) button to confirm.

##### Setting Month/Date

Initial Month/Date is 1/01, when the Month display is flashing, press the (8) button continuously, the month will increase continuously 1 by 1, press the (8) button continuously and it will reduce continuously 1 by 1, press the (1↓) button to confirm, and do in the same way to set the date, Press the (1↓) button to confirm.

##### Setting hour system

Initial hour system is 12 hour, press the (8) button or (8) button and then you can switch 12 or 24 hour system, press the (1↓) button to confirm.

##### Setting Time

When the hour display is flashing, press the (8) button continuously, the hour will increase continuously 1 by 1, press the (8) button continuously and it will reduce continuously 1 by 1, press the (1↓) button to confirm, and do in the same way to set the minute. Press the (1↓) button to confirm.

##### Setting Volume

When display with SP is flashing, press the (8) button or (8) button to switch volume 1, volume 2, volume 3 or OFF. Press the (1↓) button to confirm. After the setting is completed, the device switches off automatically and save the setting result.

### 4.3 Cuff tube connection

Insert the cuff tube into the opening on the left side of the monitor (As shown in picture-05)

## 5 Measurement Procedure


### 5.1 Before measurement

- Avoid eating and smoking as well as all forms of exertion directly before measurement. These factors influence the measurement result. Find time to relax by sitting in an armchair in a quiet atmosphere for about ten minutes before taking a measurement.
- Remove any garment that fits closely to your upper arm.
- Always measure on the same arm (normally left).

### 5.2 Fitting the Cuff

Please refer to picture-05

1. Wrap the cuff around your upper left arm. The rubber tube should be on the inside of your arm extending downward to your hand. Make certain the cuff lies approximately 2 to 3 cm above the elbow. Important! The  $\Phi$  on the edge of the cuff (Artery Mark) must lie over the artery which runs down the inner side of the arm.
2. To secure the cuff, wrap it around your arm and press the hook and loop closure together.
3. There should be little free space between your arm and the cuff. You should be able to fit 2 fingers between your arm and the cuff. Cuffs that don't fit properly result in false measurement values. Measure your arm circumference if you are not sure of proper fit.
4. Lay your arm on a table (palm upward) so the cuff is at the same height as your heart. Make sure the tube is not kinked.

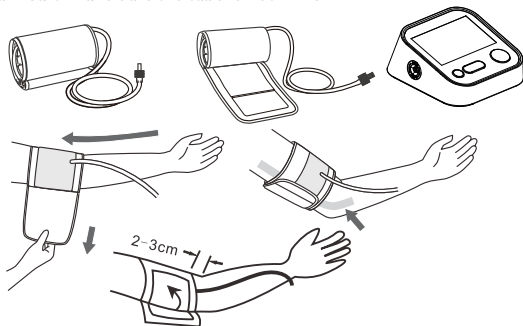

Picture-05

### 5.3 Measure Procedure

The device is designed to take measurements and store the measurement values in memory for two people using User ID 1 and User ID 2.

Refer to picture-06

1. Sit comfortably in a chair with your feet flat on the floor.
2. Select your User ID (1 or 2).
3. Stretch your arm forward on the desk and keep relaxing, make sure the palm of hand is upturned. Make sure arm is in correct position, to avoid body movement. Sit still and do not talk or move during the measurement. After the cuff has been appropriately positioned on the arm and connected to the blood pressure monitor, the measurement can begin:

#### Operate via the App on smart phone with Bluetooth

1. Install the App from Google play store or Apple app store. Open Bluetooth on smart phone, and then turn on the App.
2. Bluetooth pairing for the first time  
Press the (  $\text{Bluetooth symbol}$  ) button and the Bluetooth symbol (  $\text{Bluetooth symbol}$  ) will flash , Then operate bluetooth pairing according to the Settings on the APP, The bluetooth symbol (  $\text{Bluetooth symbol}$  ) will stop flashing after the connection is successful.
3. Bluetooth paired successfully  
Bluetooth will be automatically searched and connected when it is powered on


#### Operate on the device

1. Press the START/STOP button. The pump begins to inflate the cuff. In the display, the increasing cuff pressure is continually displayed.
2. When the device has detected your pulse, the heart symbol in the display begins to blink.
3. When the measurement has been concluded, the measured systolic and diastolic blood pressure values, as well as the pulse will be displayed.
4. The measurement results are displayed until you switch the device off. If no button is pressed for 60 seconds, the device switches off automatically.
5. Cuff self-checking symbol (  $\text{Cuff symbol}$  )  
The cuff correct symbol (  $\text{Cuff symbol}$  ) will be displayed if the cuff position is correct, otherwise the wrong symbol (  $\text{Cuff symbol}$  ) will be displayed. Please check again the cuff if the wrong symbol (  $\text{Cuff symbol}$  ) is displayed.
6. Movement error symbol (  $\text{Movement symbol}$  )  
The Movement Error Symbol (  $\text{Movement symbol}$  ) is displayed if you move your body during the measurement. Please remove the cuff, and wait 2-3 minutes. Reapply the cuff and take another measurement.

**NOTE:**

Patient Position:

1. Comfortably seated
2. Legs uncrossed
3. Feet flat on the floor
4. Back and arm supported
5. Middle of the CUFF at the level of the right atrium of the heart

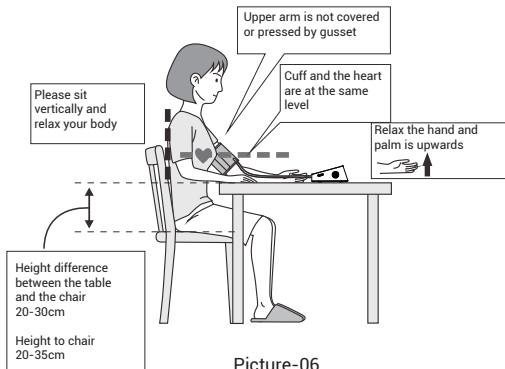

Picture-06

**5.4 Irregular Heartbeat Detector**

This symbol (♥) - indicates that certain pulse irregularities were detected during the measurement.

In this case, the result may deviate from your normal basal blood pressure – repeat the measurement.

Information for the physician on frequent appearance of the Irregular Heartbeat Symbol.

This instrument is an oscillometric blood pressure monitor device that also analyzes pulse frequency during measurement. The instrument is clinically tested.

If pulse irregularities occur during measurement, the irregular heartbeat symbol is displayed after the measurement. If the symbol appears more frequently (e.g. several times per week on measurements performed daily) or if it suddenly appears more often than usual, we recommend the patient to seek medical advice. The instrument does not replace a cardiac examination, but serves to detect pulse irregularities at an early stage.


**5.5 Error Indicates**

The following symbol will appear on the display when measuring abnormal

SYMBOL	CAUSE	CORRECTION
No display appears	Weak battery or improper placement	Replace both batteries with new ones. Check the battery installation for proper placement of the battery polarities.
Er 1	Sensor abnormal	Check if the pump is working or not. If it is working, then the problem is sensor abnormal. Please send it to the local distributor.
Er 2	Monitor could not detect pulse wave or cannot calculate the blood pressure data	Tie the cuff correctly and make the measurement again. If the error is still displayed, please send it to local distributor
Er 3	Measurement results is abnormal or out the measurable range of blood pressure	Please keep quiet and measure again
Er 4	Too loose cuff or air leakage	Tie the cuff correctly and make sure the air plug is properly inserted in the unit
Er 5	The air tube is crimped or the cuff is tied too tight.	Correct it and make the measurement again
Er 6	The sensor is sensing great fluctuation in the pressure	Please keep quiet and don't move
Er 7	The pressure that the sensor sensing is over the limit	Please send back to the local distributor
Er 8	The demarcation is incorrect or the device has not been demarcated	Please send back to the local distributor
HI	The pulse rate exceeds the upper limit (>170 per minute)	Beyond the measurement range, normal reminder
LO	The pulse rate is less than the lower limit (<40 per minute)	Beyond the measurement range, normal reminder














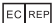




## Trouble removal

Problem	Check	Cause and solutions
No power	Check the battery power	Replace new one
	Check the polarity position	Installation for proper placement of the batteries polarities
No inflation	Whether the plug insert	Insert into the air socket tightly
	Whether the plug broken or leak	Change a new cuff
Err and stop working	Whether move the arm when inflate	Keep the body peaceful
	Check if chatting when measured	Keep quite when measure
Cuff leak	Whether the cuff wrap too loose	Wrap the cuff tightly
	Whether the cuff broken	Change a new cuff
 Please contact the distributor if you can't solve the problem, do not disassemble the unit by yourself!		

## SYMBOL DESCRIPTIONS

The following symbols may appear in this manual, on the Digital Blood Pressure Monitor F1701T, or on it's accessories. Some of the symbols represent standards and compliances associated with the Digital Blood Pressure Monitor F1701T and its use.


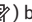
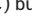
	DISPOSAL: Do not dispose this product as unsorted municipal waste. Collection of such waste separately for special treatment is necessary.
	CE Mark.
<b>IP21</b>	The degree of avoid ingress of water or particulate matter into ME equipment.
	Authorized Representative in the European Community.
	Date of manufacture.
	Manufacturer.
<b>SN</b>	Specifies serial number.
	Type BF applied part.
	Direct current.

	Follow instructions for use.
	MR unsafe.
	Medical device.
	Put up.
	Fragile.
	Afraid of the rain.
	Fear of the sun.
	Handle gently.
	Temperature range.
No Sterilize requirement.	
Not category AP/APG equipment.	
Mode of operation: continuous.	


## 5.6 Memory

Each unit stores 120 sets measurements for 2 users, totally 240 sets (User 1 and 2). Measurements for each user are stored separately. Be certain that you are viewing the measurements for the correct user.



### View the memory


With the unit off, press the () button. The monitor will display User ID and an average value of the last 3 times measurements stored in the unit. ( If measurements are less than 3 sets, directly display the first set ) Each time you press the () button, it will display the memory value from the latest to the oldest in turn. Each time you press the () button, it will display the memory value from the oldest to the latest in turn.

### Switch user memory

When viewing the memory, press the () button, it will switch the User ID, and then display another user memory value .

### Delete memory:

In average value memory viewing mode, the average value symbol (  ) is being displayed, long press the (  ) button for 3 seconds, then it will delete all measurements for the current user.

In single set memory viewing mode, long press the (  ) button for 3 seconds, then it will delete only a set measurement being displayed.

### Note:

If you decide to delete the all memory, please keep the memory in another way, incase you need it some days later. Take the battery out won't lead to a memory missing.

### 5.7 Discontinuing a Measurement

If it is necessary to interrupt a blood pressure measurement for any reason (e.g the patient feels unwell), the Start/Stop button can be pressed at any time. The device then immediately lowers the cuff pressure automatically.

### 5.8 Using the AC Adapter

You may also operate this monitor using the AC adapter (output d.c. 5V 1A with Type c connector).

1. Ensure that the AC adapter and cable are not damaged.
2. Plug the adapter cable into the AC adapter port on the right side of the blood pressure monitor.
3. Plug the adapter into your electrical outlet. When the AC adapter is connected, no battery current is consumed.

**Note:** No power is taken from the batteries while the AC adapter is connected to the monitor. If electrical power is interrupted, (e.g., by accidental removal of the AC adapter from the outlet) the monitor must be reset by removing the plug from the socket and reinserting the AC adapter connection.

## 6 Care and Maintenance

Wash hands after each time measurement.

If one device is used by different patients, wash hands before and after each use.

1. Do not expose the device to either extreme temperatures, humidity, dust or direct sunlight.
2. The cuff contains a sensitive air-tight bladder. Handle this cuff carefully and avoid all types of stress through twisting or buckling.
3. Clean the device with a soft, dry cloth. Do not use gas, thinners or similar solvents. Spots on the cuff can be removed carefully with a damp cloth and soapsuds, if necessary, 70% isopropanol can be used. The cuff with bladder must not be washed in a dishwasher, clothes washer, or submerged in water.
4. Handle the tube carefully. Do not pull on it. Do not allow the tubing to kink and keep it away from sharp edges.
5. Do not drop the monitor or treat it roughly in any way. Avoid strong vibrations.
6. Never open the monitor! This invalidates the manufacturer's warranty.
7. Batteries and electronic instruments must be disposed of in accordance with the locally applicable regulations, not with domestic waste.

### 6.1 Accuracy test

Sensitive measuring devices must be checked for accuracy from time to time. We recommend a periodical inspection of your device by an authorized dealer every 1 year. Please turn to local distributor or the manufacturer.

## 7 Warranty/Service

Your blood pressure monitor is guaranteed for 2 years against manufacturers' defects for the original purchaser only, from date of purchase. The warranty does not apply to damage caused by improper handling, accidents, unprofessional use, not following the operating instructions or alterations made to the instrument by third parties.

Warranty only applies to the main device and its cuff. All other accessories are not covered by warranty.

There are no user serviceable parts inside. Batteries or damage from old batteries is not covered by the warranty.

## 8 Certifications

Device standard:

This device is manufactured to meet the blood pressure monitors:

IEC 80601-2-30 / IEC60601-1-11 / IEC60601-1

Electromagnetic compatibility:

Device fulfills the stipulations of the International standard

IEC60601-1-2

The device was clinically investigated and the safety and efficacy is meet the requirement of ISO 81060-2. If you need to acquire a copy of the summary of the Clinical Investigation, please contact the manufacturer.

## 9 Technical Specifications

Model	F1701T
Weight	214g (batteries and AC adapter is not included)
Display	56×85mm(2.20"×3.34") LCD Digital Display
Size	121(W)×109 (L)×52(H) mm(4.76"×4.29"×2.05")
Accessories	Wide range rigid cuff 8.7" – 16.5" (22 - 42 cm)
Operating Conditions	Temperature: 5 °C to 40 °C; Humidity: 15% to 93% RH;
Storage And Shipping Conditions	Temperature: -25 °C to 70 °C; Humidity: ≤ 93% RH;
Atmospheric pressure range	70kPa~106kPa
Measuring method	Oscillometric
Pressure sensor	Resistive
Measuring range	SYS: 60-270mmHg DIA:40-220mmHg Pulse: 40 to 170 per minute
Cuff pressure display range	0-295mmHg
Memory	Automatically stores the last 120 measurements for 2 users (total 240)
Measuring resolution	1 mmHg
Accuracy	Pressure within ± 3 mmHg / pulse ± 5 % of the reading
Power source	3×AA batteries, 4.5 V AC adapter INPUT: a.c. 100-240V 50/60HZ OUTPUT: d.c. 5V 1A
Packaging list	1×Main Device, 1×Cuff, 1×Users manual,
Users	Adult
Expected service life of the device: 5 years	
Technical alterations reserved!	

## 10 EMC Declaration

The ME EQUIPMENT or ME SYSTEM is suitable for home healthcare environments and so on.

Warning: Don't near active HF surgical equipment and the RF shielded room of an ME system for magnetic resonance imaging, where the intensity of EM disturbances is high.

Warning: Use of this equipment adjacent to or stacked with other equipment should be avoided because it could result in improper operation. If such use is necessary, this equipment and the other equipment should be observed to verify that they are operating normally.

Warning: Use of accessories, transducers and cables other than those specified or provided by the manufacturer of this equipment could result in increased electromagnetic emissions or decreased electromagnetic immunity of this equipment and result in improper operation."

Warning: Portable RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 30 cm (12 inches) to any part of the Blood Pressure Monitor (F1701T), including cables specified by the manufacturer. Otherwise, degradation of the performance of this equipment could result.

If any: a list of all cables and maximum lengths of cables (if applicable), transducers and other ACCESSORIES that are replaceable by the RESPONSIBLE ORGANIZATION and that are likely to affect compliance of the ME EQUIPMENT or ME SYSTEM with the requirements of Clause 7 (EMISSIONS) and Clause 8 (IMMUNITY). ACCESSORIES may be specified either generically (e.g. shielded cable, load impedance) or specifically (e.g. by MANUFACTURER and EQUIPMENT OR TYPE REFERENCE).

If any: the performance of the ME EQUIPMENT or ME SYSTEM that was determined to be ESSENTIAL PERFORMANCE and a description of what the OPERATOR can expect if the ESSENTIAL PERFORMANCE is lost or degraded due to EM DISTURBANCES (the defined term "ESSENTIAL PERFORMANCE" need not be used).

### Technical description

1.all necessary instructions for maintaining BASIC SAFETY and ESSENTIAL PERFORMANCE with regard to electromagnetic disturbances for the excepted service life.

2. Guidance and manufacturer's declaration -electromagnetic emissions and Immunity

Table 1

Guidance and manufacturer's declaration - electromagnetic emissions	
Emissions test	Compliance
RF emissions CISPR 11	Group 1
RF emissions CISPR 11	Class B
Harmonic emissions IEC 61000-3-2	Class A
Voltage fluctuations/ flicker emissions IEC 61000-3-3	Applied

Table 2

Guidance and manufacturer's declaration - electromagnetic Immunity		
Test frequency	Modulation	IMMUNITY TEST LEVEL (A/m)
30 kHz	CW	8
134,2 kHz	Pulse modulation <sup>a)</sup> 2,1 kHz	65 <sup>c)</sup>
13,56 MHz	Pulse modulation <sup>b)</sup> 50 kHz	75 <sup>c)</sup>
<p>a) This test is applicable only to ME EQUIPMENT and ME SYSTEMS intended for use in the HOME HEALTHCARE ENVIRONMENT. b) The carrier shall be modulated using a 50 % duty cycle square wave signal. c) r.m.s., before modulation is applied.</p>		

Table 3

Guidance and manufacturer's declaration - electromagnetic Immunity		
Immunity Test	IEC 60601-1-2 Test level	Compliance level
Electrostatic discharge (ESD) IEC 61000-4-2	±8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15 kV air	±8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15 kV air
Electrical fast transient/burst IEC 61000-4-4	Power supply lines: ±2 kV input/output lines: ±1 kV 100 kHz repetition frequency	Power supply lines: ±2 kV
Surge IEC 61000-4-5	line(s) to line(s): ±0.5 kV line(s) to earth: ±2 kV line(s) to line(s): ±1 kV	line(s) to line(s): ±0.5 kV line(s) to line(s): ±1 kV.
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	0% 0.5 cycle At 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315° 0% 1 cycle And 70% 25/30 cycles Single phase: at 0 0% 300 cycle	0% 0.5 cycle At 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315° 0% 1 cycle And 70% 25/30 cycles Single phase: at 0 0% 300 cycle
Power frequency magnetic field IEC 61000-4-8	30 A/m 50Hz/60Hz	30 A/m 50Hz/60Hz
Conducted RF IEC61000-4-6	150KHz to 80MHz: 3Vrms 6Vrms (ISM and amateur radio bands) 80% Am at 1kHz	150KHz to 80MHz: 3Vrms 6Vrms (ISM and amateur radio bands) 80% Am at 1kHz
Radiated RF IEC61000-4-3	10 V/m 80 MHz – 2,7 GHz 80 % AM at 1 kHz	10 V/m 80 MHz – 2,7 GHz 80 % AM at 1 kHz
Proximity magnetic fields IEC 61000-4-39	30 kHz: 8A/m 134.2 kHz: 65A/m 13.56 MHz: 7.5A/m	30 kHz: 8A/m 134.2 kHz: 65A/m 13.56 MHz: 7.5A/m
NOTE 1: Ur is the a.c. mains voltage prior to application of the test level.		

Table 4

Guidance and manufacturer's declaration - electromagnetic Immunity					
Radiated RF IEC61000-4-3 (Test specifications for ENCLOSURE PORT IMMUNITY to RF wireless communication s equipment)	Test Frequency (MHz)	Band <sup>a)</sup> (MHz)	Service <sup>a)</sup>	Modulation <sup>a)</sup>	IMMUNITY TEST LEVEL (V/m)
	385	380-390	TETRA 400	Pulse modulation <sup>b)</sup> 18 Hz	27
	450	430-470	GMRS 460, FRS 460	FM <sup>c)</sup> ±5 kHz deviation 1 kHz sine	28
	710	704-787	LTE Band 13, 17	Pulse Modulation <sup>b)</sup> 217 Hz	9
	745				
	780				
	810				
	870	800-960	GSM 800/900, TETRA 800, iDEN 820, CDMA 850, LTE Band 5	Pulse Modulation <sup>b)</sup> 18 Hz	28
	930				
	1720	1700-1990	GSM 1800; CDMA 1900; GSM 1900; DECT; LTE Band 1, 3, 4, 25; UMTS	Pulse Modulation <sup>b)</sup> 217 Hz	28
	1845				
	1970				
	2450	2400-2570	Bluetooth, WLAN, 802.11 b/g/n, RFID 2450, LTE Band 7	Pulse Modulation <sup>b)</sup> 217 Hz	28
	5240	5100-5800	WLAN 802.11 a/n	Pulse Modulation <sup>b)</sup> 217 Hz	9
	5500				
	5785				

If necessary to achieve the IMMUNITY TEST LEVEL, the distance between the transmitting antenna and the ME EQUIPMENT or ME SYSTEM may be reduced to 1 m. The 1 m test distance is permitted by IEC 61000-4-3.

a) For some services, only the uplink frequencies are included.  
b) The carrier shall be modulated using a 50 % duty cycle square wave signal.  
c) As an alternative to FM modulation, the carrier may be pulse modulated using a 50 % duty cycle square wave signal at 18 Hz. While it does not represent actual modulation, it would be worst case.